

Fabrifoam® introduces a revolutionary material with greater support.

# X-TREME™




When semi-rigid support is the protocol, Fabrifoam® X-treme™ is the choice.

## FABRIFOAM® X-TREME™

Bracing, compression and cushioning for a new level of comfort and excellence in musculoskeletal care.


- Advanced triple composite version of Fabrifoam provides the ultimate in moderate to semi-rigid support and compression
- Breathable, non-migratory, latex-free & neoprene-free
- Moisture-wicking, hypoallergenic, washable, trimmable
- Easy on and off with hook & loop closures

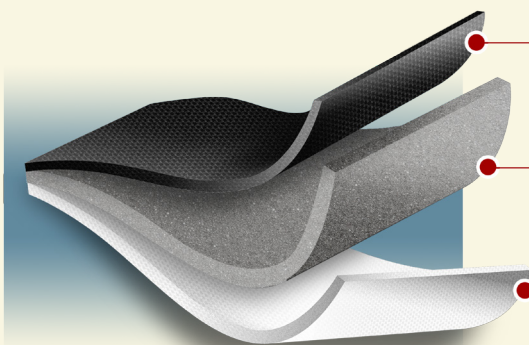

- 1 **Fabric Layer** A brushed, nylon/spandex knitted blend with various stretch capabilities. Provides adjustability for optimal fit and ease of use.
- 2 **X-layer Core** Ester-based, high elasticity polyurethane foam with durable low compression-set properties for enhanced stability and support.
- 3 **Foam Layer** Polyether-polyurethane open-cell foam with a high coefficient of friction for slip resistance, breathability and comfortable 24/7 skin contact.

## ORIGINAL FABRIFOAM®

Fabrifoam® is a fabric/foam composite marrying an open-cell, hydrophilic foam with high quality medical grade fabrics specifically engineered and designed for the treatment and prevention of musculoskeletal conditions.

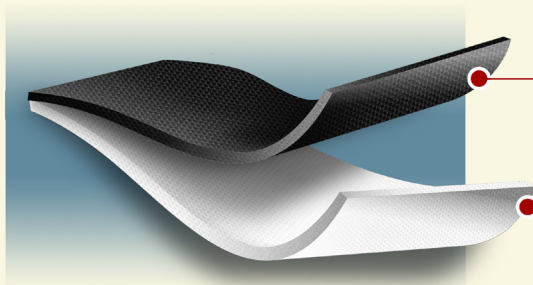

- 1 **Fabric Layer** A brushed, nylon/spandex knitted blend with various stretch capabilities. Provides adjustability for optimal fit and ease of use.
- 2 **Foam Layer** Polyether-polyurethane open-cell foam with a high coefficient of friction for slip resistance, breathability and comfortable 24/7 skin contact.



# FABRIFOAM® X-TREME™ PRODUCTS

Stronger. Greater Compression. Better Cushioning. No Uncomfortable Bulk.



## Thumb Spica X-treme™

For CMC and MCP joint instability and injury. Provides flexible support to aid in restricting movement of the thumb. Helps decrease pain and swelling.



## CarpalGard™ X-treme™

Ideal for Carpal Tunnel Syndrome, wrist pain and strain and other wrist ailments.



## ElbowGard™ X-treme™

Recommended for tennis elbow, golfer's arm and tendonitis.



## AnkleGard™ X-treme™

Strong, yet flexible and ideal for ankle sprains, pain, injury and instability. Helpful for post-op recovery.



## KneeGard™ X-treme™

Helps with knee sprains, runner's knee, injury prevention, and rehab.



## PSC™ X-treme™

Dynamic orthosis to treat plantar fasciitis, excessive pronation, chronic heel pain, heel spur syndrome and shin splints. Reduces the heel strike force and biomechanically moves the foot into a more neutral alignment. Simulates low-dye taping, simple to use, with or without orthotics."



## PattStrap™ X-treme™

Wear above the knee to relieve the stress and pain of Iliotibial Band Syndrome. Wear below the knee to treat patellar tendonitis, Osgood Schlatter Disease and chondromalacia.

